


Science made smarter

AS608 / AS608e

Quick hearing
screening
made real



Diagnostic Suite software
(AS608e only)

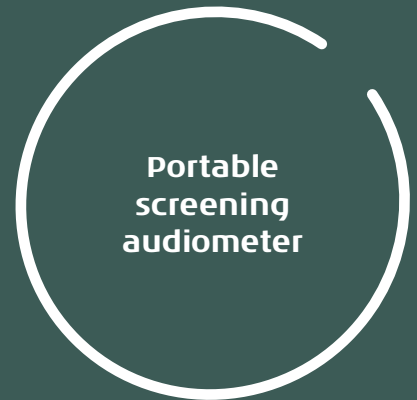

Portable
screening
audiometer



Interacoustics

Audiometry
~~~~~

Tympanometry

ABR

OAE

Hearing Aid Fitting

Balance



# Portable and easy to use

**The AS608 / AS608e portable audiometer is ideal for schools, industry and primary medical practices or wherever quick hearing screening evaluations are performed. Easy to use, it requires minimal training prior to use and offers calibrated pure tone and warble tone stimuli over a wide range of frequencies.**

#### **Portable**

The AS608 / AS608e includes a dedicated, lightweight carrying bag with shoulder strap that will accommodate the AS608 / AS608e, headset, and audiogram charts.

#### **Batteries or power supply**

The AS608 / AS608e will operate with 3 standard AA batteries or a medically approved external power supply. Typical battery life is six months.

#### **PC-integration (AS608e only)**

The AS608e offers PC-integration through a USB cable using the Diagnostic Suite Windows® software for transferring audiometric data to the PC. The diagnostic suite offers features such as:

- Print Wizard
- Advanced report page
- Session handling
- Database storage: OtoAccess™ and Noah compatible

#### **Auto Threshold (AS608e only)**

The AS608e includes the Hughson Westlake automatic pure tone test procedure. When the test is completed the results are easily recalled from the internal memory of the AS608e and displayed in the Diagnostic Suite PC software.

#### **A preview of the benefits**

- Robust construction
- 11 frequencies from 125Hz to 8kHz
- Customize frequency test selections
- Output levels from -10dB to 100dBHL
- Pure tone, warble tone and pulsed tone
- Silent stimulus switch
- Battery status indicator
- Carrying bag
- Patient response switch (included with AS608e)
- External power supply (optional)
- Noise reducing headset (optional)
- Talk Forward (only AS608e)
- PC-integration (only AS608e)
- Hughson Westlake automatic pure tone test procedure (only AS608e)


Total weight  
less than  
1.6 kg/3.5 lbs



The soft carrying case



Report page from the Diagnostic Suite (AS608e only)

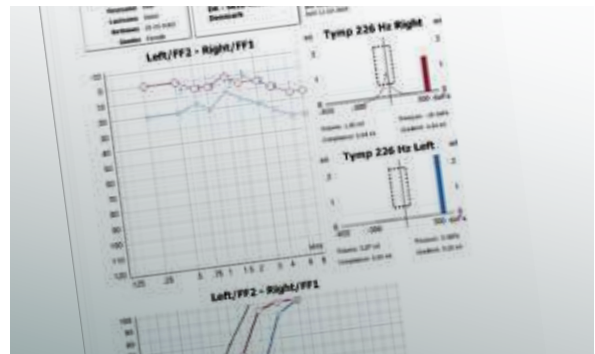

Print Layout from the Diagnostic Suite (AS608e only)

# Science made smarter

## Interacoustics is more than state-of-the-art solutions

Our mission is clear. We want to lead the way in audiology and balance by translating complexity into clarity:

- Challenges made into clear solutions
- Knowledge made practical
- Invisible medical conditions made tangible and treatable

Our advanced technology and sophisticated solutions ease the lives of healthcare professionals.

We will continue to set the standard for an entire industry. Not for the sake of science. But for the sake of enabling professionals to provide excellent treatment for their millions of patients across the globe.

[Interacoustics.com](http://Interacoustics.com)

Interacoustics A/S

Audiometer Allé 1  
5500 Middelfart  
Denmark

+45 6371 3555  
[info@interacoustics.com](mailto:info@interacoustics.com)

[interacoustics.com](http://interacoustics.com)

Go online to  
explore our  
full product  
range

## Related products



**MT10**  
Handheld tympanometer



**AD226**  
Diagnostic audiometer



**Viot™**  
Video otoscope

### Product specifications

All technical and hardware specifications concerning all products can be downloaded from our website.



**Interacoustics**